


# FilmArray® Meningitis/Encephalitis Panel


1 Test. 14 Targets. All in about an hour.



Bacteria

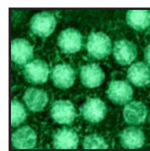

Viruses

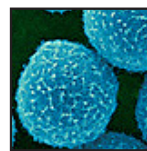


Fungi

*Escherichia coli* K1  
*Haemophilus influenzae*  
*Listeria monocytogenes*  
*Neisseria meningitidis*  
*Streptococcus agalactiae*  
*Streptococcus pneumoniae*

Cytomegalovirus (CMV)  
Enterovirus  
Herpes simplex virus 1 (HSV-1)  
Herpes simplex virus 2 (HSV-2)  
Human herpesvirus 6 (HHV-6)  
Human parechovirus  
Varicella zoster virus (VZV)

*Cryptococcus neoformans/gattii*

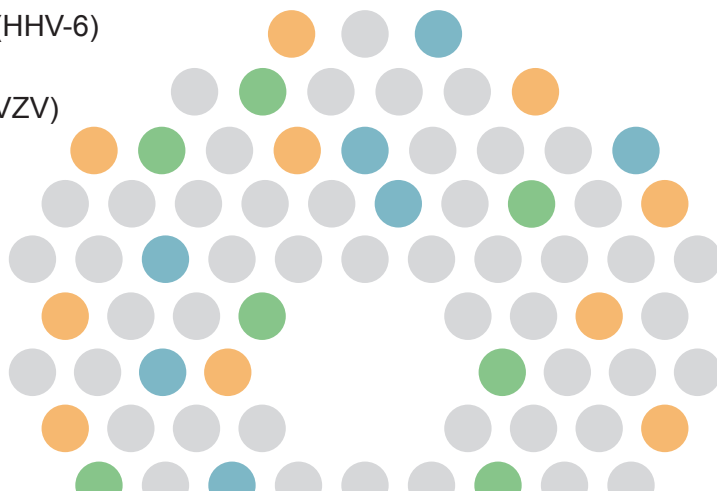

14

Targets

## Comprehensive Panel of 14 Targets

The FilmArray Meningitis/Encephalitis (ME) Panel tests cerebrospinal fluid (CSF) for a variety of pathogens including bacteria, viruses, and fungi. The integrated FilmArray system yields sample to results in about an hour, with only 2 minutes of hands-on time.

- **Simple:** 2 minutes of hands-on time
- **Easy:** No precise measuring or pipetting required
- **Fast:** Turnaround time of about 1 hour
- **Comprehensive:** 14 bacterial, viral, and fungal targets

For In Vitro Diagnostic Use  
FDA cleared and CE marked



If you are interested in a free, no obligation demonstration of the FilmArray in your laboratory visit [www.BioFireDX.com](http://www.BioFireDX.com) or call 1-800-735-6544.

FREE  
Demo!



FilmArray®

# Panel Specifications

## Sample Requirements

- Sample Type: Cerebrospinal Fluid (CSF)
- Sample Volume: 200 µL

## Performance Parameters

- Hands-on time: Approx. 2 minutes
- Run turnaround time: About 1 hour


## How Does the FilmArray Work?

The FilmArray reagent pouch stores all the necessary reagents for sample preparation, reverse transcription, PCR and detection in a freeze-dried format. Prior to a run, the user injects hydration solution and sample combined with sample buffer into the pouch. The FilmArray instrument does the rest.


First, the FilmArray extracts and purifies all nucleic acids from the sample. Next, the FilmArray performs a nested multiplex PCR. During the first-stage PCR, the FilmArray performs a single, large volume, massively multiplexed reaction. Last, individual singleplex second-stage PCR reactions detect the products from the first stage PCR.

Using endpoint melting curve data, the FilmArray software automatically generates a result for each target in a single report.

## Setting up the FilmArray is Easy – Sample in, Results out

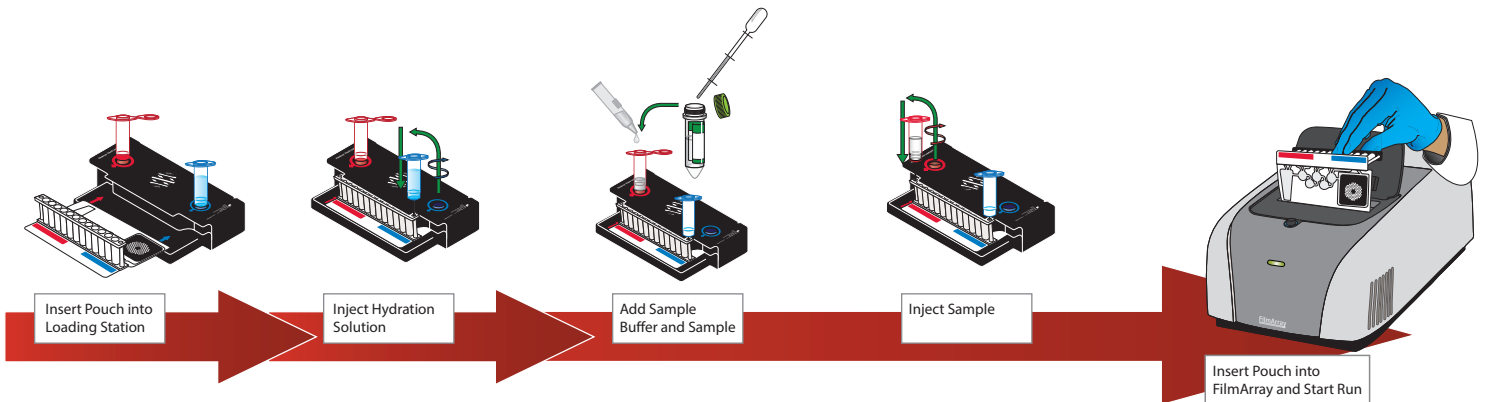

## The FilmArray Pouch and Analysis Report

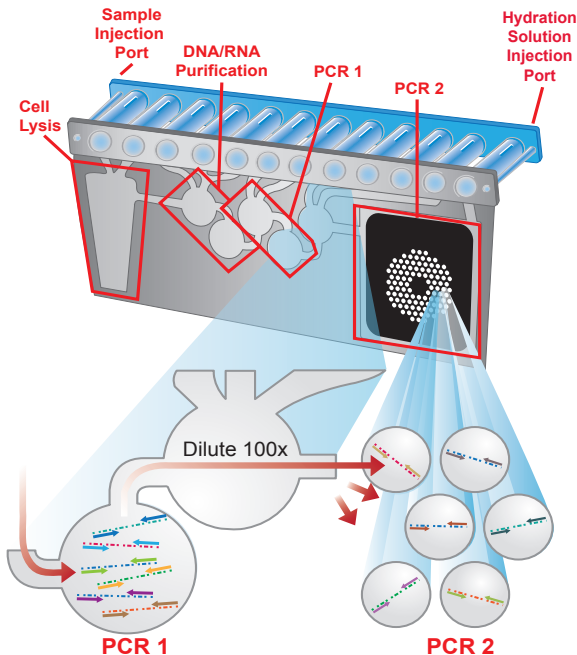

FilmArray® Meningitis / Encephalitis (ME) Panel - IVD		BIO FIRE	
<a href="http://www.BioFireDx.com">www.BioFireDx.com</a>			
<b>Run Summary</b>			
Sample ID:	14067-0295	Run Date:	21 Jun 2015 1:47 PM
Detected:	Human herpesvirus 6	Controls:	Passed
<p><b>WARNING:</b> The FilmArray ME panel does not distinguish between latent and active CMV and HHV-6 infections. Detection of these viruses may indicate primary infection, secondary reactivation, or the presence of latent virus. Results should always be interpreted in conjunction with other clinical, laboratory and epidemiological data.</p>			
<b>Result Summary</b>			
		<b>Bacteria</b>	
Not Detected	<i>Escherichia coli</i> K1		
Not Detected	<i>Haemophilus influenzae</i>		
Not Detected	<i>Listeria monocytogenes</i>		
Not Detected	<i>Neisseria meningitidis</i>		
Not Detected	<i>Streptococcus agalactiae</i>		
Not Detected	<i>Streptococcus pneumoniae</i>		
		<b>Viruses</b>	
Not Detected	Cytomegalovirus		
Not Detected	Enterovirus		
Not Detected	Herpes simplex virus 1		
Not Detected	Herpes simplex virus 2		
Not Detected	Human herpesvirus 6		
Not Detected	Human parechovirus		
✓ Detected	Varicella zoster virus		
		<b>Yeast</b>	
Not Detected	<i>Cryptococcus neoformans/gattii</i>		
<b>Run Details</b>			
Pouch:	ME Panel v1.4	Protocol:	CSF v2.0
Run Status:	Completed	Operator:	Ashley Hunter (ah)
Serial No.:	01422247	Instrument:	ITI FA "FA2430"
Lot No.:	140411A		

The purchase of FilmArray System includes a limited, non-transferable license under U.S. Patent No. 5,871,908, owned by Evotec Biosystems GmbH and licensed to Roche Diagnostics GmbH, to use only the enclosed amount of product according to the specified protocols. No right is conveyed, expressly, by implication, or by estoppel, to use any instrument or system under any claim of U.S. Patent No. 5,871,908, other than for the amount of product contained herein.

The information contained in this document is subject to change without notice. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of BioFire Diagnostics, LLC. FilmArray

Software, Detector, and Metacall software modules ©2002–2015 BioFire Diagnostics, LLC. BioFire Diagnostics, BioFire, the BioFire logo, FilmArray and LCGreen are trademarks of BioFire Diagnostics, LLC or BioFire Defense, LLC and are registered trademarks of BioFire Diagnostics, LLC or BioFire Defense, LLC in the United States. All other names of products and brands appearing in this manual are trademarks or registered trademarks of their respective owners.

The purchase of this product includes a limited, nontransferable license under specific claims of one or more U.S. patents as listed on BioFire Diagnostics' Web site (<http://www.biofiredx.com/LegalNotices/>) and owned by BioFire and the University of Utah Research Foundation.